


ELECTRODES FOR EMISSION SPECTROGRAPHY

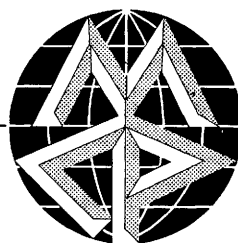

High purity

Rods — in 6 grades

Preformed electrodes

Silver and iron comparison electrodes

Discs · Powder

Catalogue on request

***** *IMMEDIATE DELIVERY*

from

Mining & Chemical Products Limited

86 STRAND, LONDON, W.C.2

TEL: COVENT GARDEN 3393


AND

70-76 ALCESTER ROAD SOUTH, KINGS HEATH

BIRMINGHAM 14

TEL: Highbury 2281

Pyrex electrically conducting coated tubing

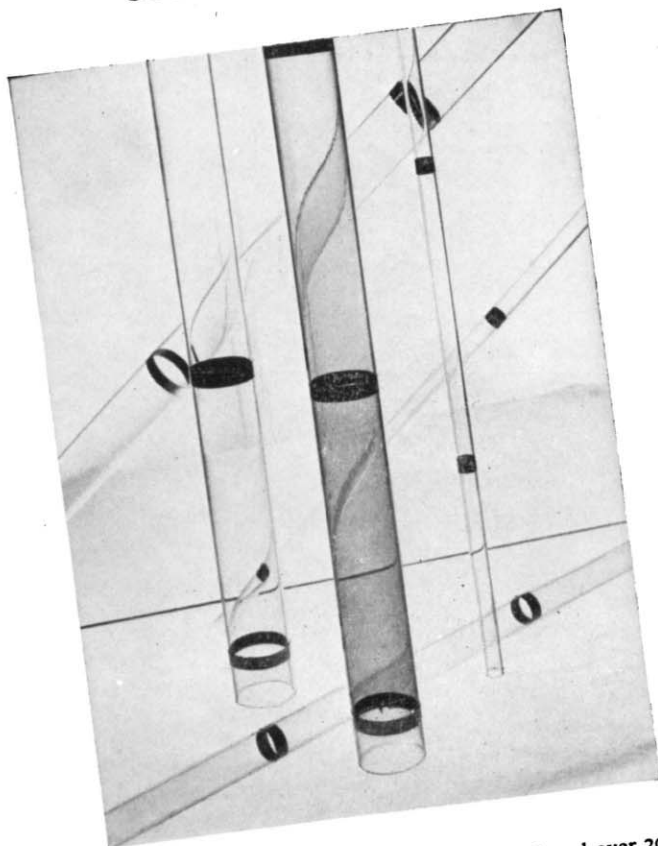

Pyrex E.C. (Electrically Conducting) glass was first developed over 20 years ago in the Research Laboratories of the Corning Glass Works, Corning, New York, U.S.A. By combining the electrical conductivity of metals with the transparency and chemical stability of glass, a new approach to process heating was possible.

After full investigation by our Research and Development staff, this product has been manufactured and made available in this country by James A. Jobling & Co. Limited.

Description

E.C. coated tubes are made from the well-known Pyrex borosilicate heat-resisting glass. A microscopically thin transparent film of metal is permanently bonded to the outer surface of a glass tube to form an electrically conducting skin.

The tubes are manufactured in standard sizes of 19, 30 and 64mm. outside diameter and 30 and 48 inch lengths in each size. To facilitate the control of the

read the full story . . .

in this month's **'PYREX GLASS'**

This is the 'Article of the Month' in the current issue of 'Pyrex Glass', a technical bulletin distributed free by James A. Jobling to all who are interested in the advances now being made in the development of Laboratory and Scientific glassware. 'Pyrex Glass' is published every two months, and each issue includes specially written articles.

TO RECEIVE FREE COPIES OF 'PYREX GLASS' REGULARLY

please complete this coupon and post it to us

'Pyrex' is the registered trade mark of James A. Jobling & Co. Ltd.


To The Industrial Sales Department, James A. Jobling & Co. Ltd.,
Wear Glass Works, Sunderland.

Please send copies of this issue of 'Pyrex Glass' and of all
future issues, to:

name (and department).....

address


.....

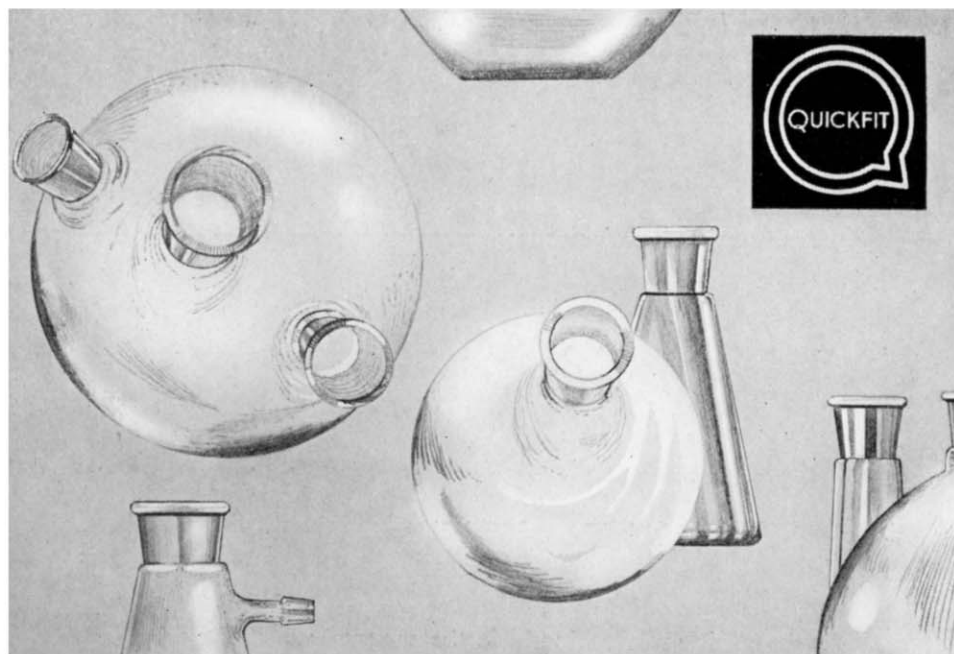

The Comprehensiveness of "QUICKFIT" leads to progressive economy in laboratory glassware

Comprehensiveness is the keynote of "Quickfit" interchangeable laboratory glassware. Your flask requirements must be very unusual if you cannot find them within the "Quickfit" range. The same may be said of adapters, stoppers and extraction apparatus, as well as all

the other types of laboratory glassware made by "Quickfit." That name is the synonym for comprehensive laboratory glassware. It brings complete standardisation to laboratory glass equipment. It enables you to do more with rather less. "Quickfit" is the basis of progressive laboratory economy.

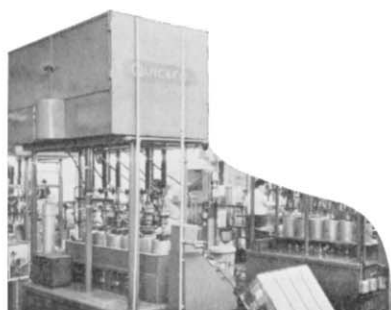

Semi-automatic machines grinding "Quickfit" precision ground joints

QUICKFIT & QUARTZ Ltd.

Dept. O(h), "Quickfit" Works,
Heart of STONE, Staffordshire
Phone: STONE 481

Water Analysis



SP.600 Spectrophotometer


To all laboratories whose concern is the quality control of a product where the use of water is a vital part of the process — in Brewing & Distilling, Dyeing, the Manufacture of Foodstuffs, or where the impurities in boiler feed water or effluent must be accurately assessed — the analysis of the many organic and inorganic compounds in water is essential.

The Unicam SP.600 Spectrophotometer is a precise economical means of maintaining a rigid standard and of avoiding needless waste. The SP.600 will give accurate information quickly; it is simple to operate and its cost is comparatively low.


UNICAM SP.600

For full information on the SP.600 please write to

UNICAM INSTRUMENTS LTD • ARBURY WORKS • CAMBRIDGE

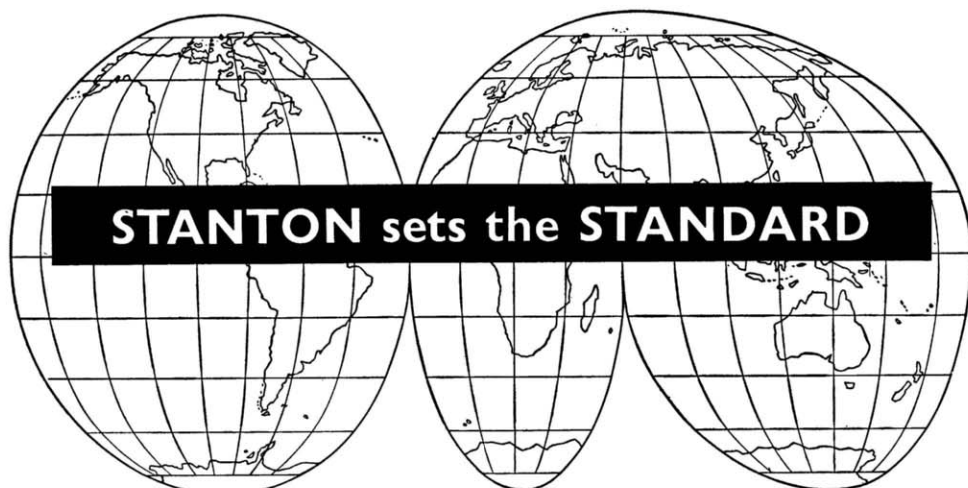

1 SECOND IN 6 YEARS corresponds to a Standard Deviation of 1 microgramme on a 200g balance. Such is the remarkable accuracy attained on a Stanton High Precision balance. That's why the laboratories in 5 continents, charged with maintaining the world standards of mass, turn to Stanton for their precision balances and weights.

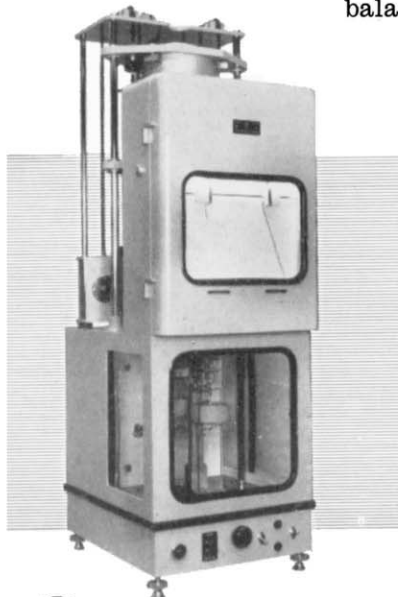

That is why the

STANTON

THERMO RECORDING BALANCE
has achieved its success

The THERMO-BALANCE automatically records weight, time and temperature changes. Air-damped, suitable for total load of 50g; may be supplied with sensitivity of 1mg or 1/10mg. Standard furnace (illustrated) suitable for 1000°C, but alternative furnaces for both high and low temperatures can be made. Once set, the Thermo-Balance will reproduce and record weight changes, time and temperature for periods up to several days.

FULLY DESCRIPTIVE LEAFLET AVAILABLE



Regd.
Trade Mark

STANTON *Thermo-recording Balances*

STANTON INSTRUMENTS LTD · 119 OXFORD STREET · LONDON W.1

Telephone: Gerrard 7533 · Cables: Stanbal London



Atomic-absorption Spectroscopy

Atomic-absorption spectroscopy, originally developed by Dr. A. Walsh of the C.S.I.R.O., Melbourne, Australia, has certainly made its mark on the literature. We give here, as a matter of interest, some of the bibliography on the subject. We regret that we are unable to give a complete bibliography in this space or to supply reprints of these papers.

ANALYST:

The Determination of Exchangeable Sodium, Potassium, Calcium and Magnesium in Soils by Atomic-Absorption Spectrophotometry. *David, D. J.* **85**, 495 (1960)

The Application of Atomic Absorption to Chemical Analysis. A Review. *David, D. J.* **85**, 779 (1960)

Atomic-Absorption Spectrophotometry with Special Reference to the Determination of Magnesium. *Allan, J. E.* **83**, 466 (1958)

Determination of Zinc and Other Elements in Plants by Atomic-Absorption Spectroscopy. *David, J. E.* **83**, 655 (1958)

The Quantitative Determination of Some Noble Metals by Atomic-Absorption Spectroscopy. *Lockyer, R., Hames, G. E.* **84**, 385 (1959)

Determination of Calcium in Plant Material by Atomic-Absorption Spectrophotometry. *David, D. J.* **84**, 536 (1959)

Determination of Zinc in Metallurgical Materials by Atomic-Absorption Spectroscopy. *Gidley, J. A. F., Jones, J. T.* **85**, 249 (1960)

NATURE:

Atomic-Absorption Spectrophotometric Determination of Molybdenum and Strontium. *David, D. J.* **187**, 1109 (1960)

Determination of Nickel and Cobalt by Atomic-Absorption. *Allen, J. E.* **187**, 1110 (1960)

Determination of Silver in Lead Concentrates by Atomic-Absorption Spectroscopy. *Rawling, B. S. et al.* **188**, 137 (1960)

Determination of Magnesium in Blood Serum by Atomic-Absorption Spectroscopy. *Willis, J. B.* **184** (468I), **187** (1959)

Some Atomic Reactions by Absorption Spectroscopy. *Broida, H. P., Schiff, H. I., Sugden, T. M.* **185**, 759 (1960)

Determination of Calcium in Blood Serum by Atomic-Absorption Spectroscopy. *Willis, J. H.* **186** (4720), 249 (1960)

ANALYTICAL CHEMISTRY:

Atomic-Absorption Spectroscopy. (Survey). *Robinson, J. W.* **32**, 17A (1960)

A Study of Atomic-Absorption Spectroscopy. *Menzies, A. C.* **32**, 898 (1960)


*For details of Hilger and Watts
atomic-absorption equipment, please write for catalogue CH407/A9*

HILGER AND WATTS LTD · 98 ST PANCRAS WAY · LONDON NW1



DERRITRON *Electronic* GROUP

BEME TELECOMMUNICATIONS LTD · BEULAH ELECTRONICS LTD
 CHAPMAN ULTRASONICS LTD · DORAN INSTRUMENT CO. LTD
 ELECTRONIC VIBRATORS LTD · L.S.B. COMPONENTS LTD · RESLOSOUND LTD

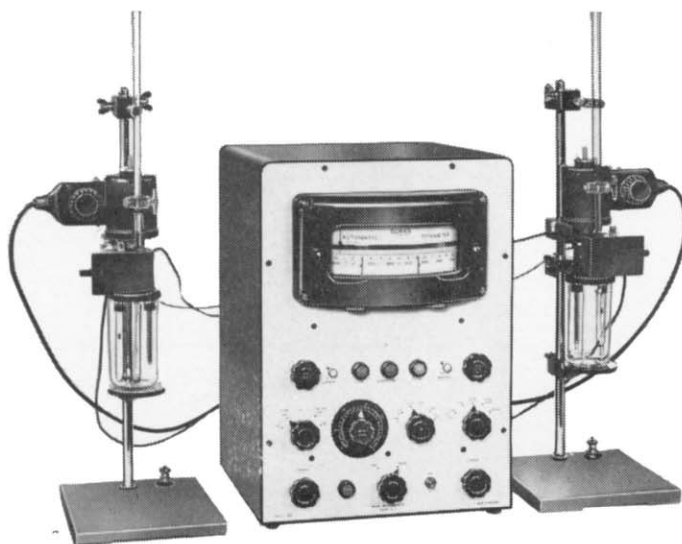

AUTOMATIC TITRIMETER

This is a continuously indicating instrument using a single staged balanced amplifier ensuring simplicity of operation. The titration may be automatic or manual and the progress of the titration can be observed visually. For further details see List 121.

SPECIFICATION

Cat. No:	M.4910
Ranges	0-14 pH and 200-0-1200 mVs.
Stability of End Point Control	Better than 0.1 pH or 10 mV.
Range of anticipation	12 pH or 1200 mVs.
Titration accuracy	± 1 drop
Rate of Reagent Flow	Max. 0.5 ml. per sec. Min. 0.01 ml. per sec.

DORAN INSTRUMENT CO. LTD

Please write today for fully illustrated details to:—

24 UPPER BROOK ST., MAYFAIR, LONDON, W.1. Phone: HYDE PARK 2291

A REMARKABLE INSTRUMENT

New standard of performance in colorimetric analysis

The accumulated experience gained during the last twenty-one years in producing photo-electric apparatus for colorimetric analysis has gone into the design and construction of the new "EEL" Spectra.

This is an instrument with all the advantages of versatility and robustness, combined with accuracy and simplicity of operation.

The main features in the unique construction of "EEL" Spectra are the high quality continuous spectrum wedge having a narrow band width giving accurate wave length selection over the range 400-700 microns and a sliding carriage to accommodate up to four cells or test tubes.

The cells can be of the order of 2.5, 5, 10 or 20 mm., the test tubes $\frac{1}{2}$ in.

The light from a prefocus lamp is passed through a heat-absorbing filter and condenser lens assembly to focus on a slit in front of the continuous spectrum wedge. A cylindrical lens is interposed between the wedge and the sample and restricts the vertical spread of the beam. Transmitted light from the sample is collected by a single standard "EEL" barrier layer photocell which is coupled to a built-in, sealed taut suspension galvanometer fitted with zero and clamping controls.

Sample solutions contained in the fused glass cells or available test tubes are mounted on a sliding carriage moving between stops. Provision is made for the containers holding

the solution under test and the reference blank to be inserted simultaneously and these may be positioned in the light-beam as required. Immediate readings are given on the galvanometer scale which is clearly marked in inverse logarithmic divisions.

The "EEL" Spectra is attractively and functionally designed with many outstanding features. Amongst the diversity of applications for which it is suited in a wide field of industrial analysis we include the following as a representative selection.

AGRICULTURE. Measurement of trace elements in soils, plants, etc.

FOOD MANUFACTURE AND BREWING. Trace elements estimation. Liquid colour analysis and control.

MEDICINE. Biochemical analysis.


METAL AND CERAMICS. Analysis of metal alloys, silicate materials, trace determinations, etc.

WATER AND SEWAGE TREATMENT. Trace determinations and chloride control.

INDUSTRIAL EDUCATION. Wide range of colorimetric tests.

Further details and price of this remarkable instrument are available from the manufacturers, EVANS ELECTROSELENIUM LTD., Sales Division, 12, St. Andrews Works, Halstead, Essex, who will be pleased to arrange a free demonstration of the "EEL" Spectra at your laboratory.

Sales and Servicing agents throughout the world.

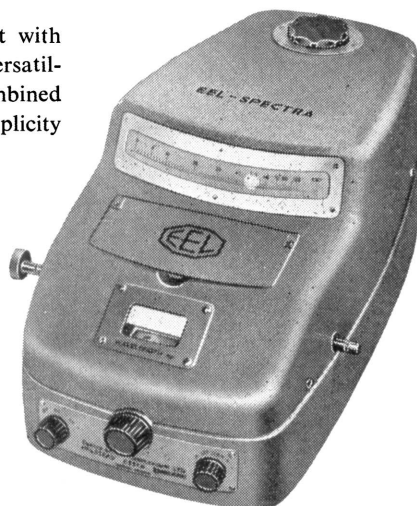

THE "EEL" SPECTRA

Choose your Cellulosic Ion-Exchangers from these Whatman trial packets

Pack 1	CELLULOSIC ION-EXCHANGE PAPERS 2 sheets 9" x 11" of each grade
Pack 2	CELLULOSIC ION-EXCHANGE POWDERS 15 g of each grade
Pack 3	CELLULOSIC ION-EXCHANGE FLOCS 15 g of each grade

The technique of ion-exchange using modified cellulose materials is attracting many new applications, many new users. To simplify your task of choosing which of the grades of Whatman floccs, powders or papers is best suited to your immediate needs, three trial packets of materials have been prepared. Each contains a good selection of the various cation and anion exchangers including cellulose phosphate, carboxylated cellulose, aminoethyl cellulose, DEAE-cellulose and Ectoelacellulose. Each is available at a reasonable price. A General Information leaflet describing the range of Whatman modified cellulose ion-exchangers is also available from

H. REEVE ANGEL & CO LTD
Gt. Britain 9 Bridewell Place London EC4
U.S.A. 9 Bridewell Place Clifton New Jersey

Sole Agents for

WHATMAN MATERIALS
Manufactured by W & R Balston Ltd

But my dear chap,

all the best people

use Oertling

nowadays



For speed, precision

and long life, Oertling

riderless micro-balances

are to-day the world standard

by which others are judged. They are backed by

a maintenance service which is also the finest

available. In the complete range,

sensitivities run from 0.02 mg.

to 0.002 mg. per division,

direct reading.



Have you sent for the latest catalogue ?

OERTLING MICRO-BALANCES

L. Oertling Ltd., Cray Valley Works, Orpington, Kent

Telephone : Orpington 25771

Overseas agents in Australia, Brazil, Canada, Holland, India, Pakistan, Union of S. Africa and U.S.A. provide sales and service. Enquiries should be sent in the first place to Head Office at the address above.

Looking for traces?



MERVYN-HARWELL SQUARE WAVE POLAROGRAPH

Mk. III

**A basic scientific advance in the field of Analysis,
Testing and Quality control.**

Whether you are "Looking for Traces" in processing Metals or Foodstuffs, carrying out research on Gas or Petroleum, or linking quality-testing with production processes, you should invest in the MERVYN-HARWELL SQUARE WAVE POLAROGRAPH

- Analyses minute quantities.
- Extremely sensitive.
- Exceptionally selective.
- Completely self-contained.
- Versatile, Adaptable.
- Operators need not be specialists.
- Strikingly convenient, mobile.
- Controls located and ganged for rapid, accurate operation.
- Built-in cell, thermostat.
- Deoxygenated supply and recorder.
- Reversible and irreversible reactions.
- Resolves 20 mv.

*Write for further information
mentioning "Analyst" to:*

MERVYN INSTRUMENTS

WOKING · SURREY · ENGLAND

WOKING 5200